

A-dec 2000



 a-dec[®]

The A-dec Difference.

The best way to discover what dentists need in a chair and delivery system is to ask. We did. And the A-dec 200 design reflects the way dental teams work, and purposefully supports efficiency in your practice. Everything about the A-dec 200 reflects the skill and quality of your care.



Backed by a 5-year warranty.





A white dental sink with a control panel and various dental handpieces. The control panel is located on the right side of the sink and features a digital display and several buttons. The handpieces are mounted on a tray below the sink. The text "Designed for Durability." is overlaid on the sink.

Designed for Durability.



Pure light. Pure dependability. “Good enough” just isn’t. A-dec designs and tests our lights for a minimum 20-year life.

Your equipment needs to work when you do. That’s why A-dec manages the manufacturing processes, so we can provide doctors with the most durable, reliable equipment. For example, we build and test each chair cylinder from raw material to the finished assembly, and back it with a 10-year warranty. And we manufacture our own high-quality delivery system tubing. Our production systems are consistently monitored to ensure that they meet A-dec quality standards. With the A-dec 200 solution, you get a strong, robust system designed to withstand years of constant use.





Designed to Perform.

A-dec 200 is a system that's not only made to last, but made to grow and evolve with your practice. Detailed features such as a stain-resistant cuspidor that rotates toward the patient, a telescoping assistant's arm for close work, and a touchpad that integrates chair, light, and cuspidor functions all add up to a better workflow. For future expansion, an open platform leaves room to add or change ancillary devices. You get peak performance now, and years from now.



- Three axes light head lets you position any A-dec operator light at just the right angle for you and your patient.
- Four preset chair positions give you added flexibility for positioning the patient.
- Ancillary accommodation allows you to add new technology as your practice evolves.
- Two-position armrests and seat articulation provide patient comfort.
- Contoured covers and minimal seams support quick and easy cleaning.

Delivery System

Keep everything you need within reach to increase the speed of treatment.



A-dec 200



A-dec 300 Continental



A-dec 300 Traditional

Touchpad

Fingertip access simplifies chair adjustments and provides quick access to dental light and cuspidor control. Optional deluxe touchpad adds controls for integrated clinical devices such as electric motors, cameras, and scalers.



Standard



Deluxe

Foot Controls

Disc style with chip blower—or new A-dec lever style for precise modulation of electric and pneumatic handpieces.



Disc foot control with chip blower



Lever foot control

Dental Lights

A-dec dental lights provide true light for accurate diagnosis, reduced eyestrain and unparalleled ergonomics.



A-dec 200 Halogen



A-dec 300 LED



A-dec 500 LED

Monitor Mount

The monitor mount positions easily over the patient for convenient patient consultation and education.



Monitor Mount

Cuspidor and Support Center

The porcelain cuspidor swivels for just the right patient position. Easy-to-fill, two-liter water bottle provides a continuous water supply through multiple procedures. The support center offers ample storage room to add ancillaries such as amalgam separators and vacuum systems.



Cuspidor and support center

Assistant's Instrumentation

Three-position holder and telescoping arm conveniently positions vacuum instruments for accessibility.



A-dec (aluminum) with touchpad



Durr (polymer) without touchpad

Formed Upholstery

A sleek vacuum-formed fit with no visible seams make cleaning and maintenance easy. See your authorized A-dec dealer for accurate color samples and the most current product information.




sky blue



indigo



sapphire

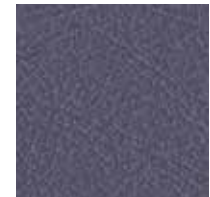


amethyst



lemongrass



black

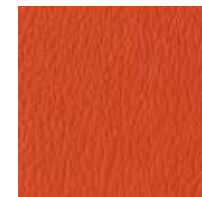

paprika



tawny

Specifications

A-dec 200 Dental Chair

Range of lift height	15.75" to 28.75" (400 mm to 730 mm)
Flexible backrest	2.25" (55 mm), seamless
Drive system	Hydraulic base and tilt
Chair back and base limit switch	Yes, stops downward movement when obstruction is detected. Patient comfort toeboard articulation
Headrest	Dual-articulating (knob style)
Armrest	Pivoting
Integrated 300-watt power supply	Yes
Foot switch	Optional
Mounting capability	Floating - 2" post mount
Upholstery options	Formed

Delivery Systems

	A-dec 200	A-dec 300 Continental	A-dec 300 Traditional
Control block	Standard 4-position	Standard 3-position, optional 4-position	Standard 3-position, optional 4-position
Left/right conversion	Yes	Yes	Yes
Tubing options	Vinyl or silicone	Vinyl or silicone	Vinyl or silicone
Handpiece holder positions	4	4	4
Brake handle(s)	N/A	1 standard, 2nd optional	1 standard, 2nd optional
Touchpad options	None, standard or deluxe	None, standard or deluxe	None, standard or deluxe
Quad-voltage intraoral light source	Optional	Optional	Optional
Foot control	Disc with chip blower or lever	Disc, disc with chip blower or lever	Disc, disc with chip blower or lever
Arm system	Balanced flexarm	Balanced flexarm with air brake A-dec 335 - support mount	Balanced flexarm with air brake A-dec 334 - support mount





Building Relationships.

A-dec is headquartered in Newberg, Oregon USA, but our customers span six continents—every one except Antarctica. Schools and institutions in more than 85 countries rely on A-dec equipment to withstand the rigors of continual use day after day, year after year. Yet with A-dec, you get more than quality equipment. You get full service and support long after the sale and installation. You get local dealers and service professionals, who will respond to your needs and keep you working. And, you get the benefit of more than 50 years of experience to help your practice grow. If you ever have a question, concern or even an idea for improvement, we're just a phone call away. We're here to help you succeed, because at A-dec we don't just build equipment. We build relationships.



A-dec Headquarters

2601 Crestview Drive
Newberg, Oregon 97132 USA
Tel: 1.800.547.1883 within USA/CAN
Tel: +1.503.538.7478 outside USA/CAN
Fax: +1.503.538.0276
a-dec.com

A-dec Australia

Tel: 1.800.225.010 within AUS
Tel: +61.(0)2.8332.4000 outside AUS

A-dec China

Tel: 400.600.5434 within CN
Tel: +86.571.89026088 outside CN

A-dec United Kingdom

Tel: 0800.233285 within UK
Tel: +44.(0).24.7635.0901 outside UK

A-dec, the A-dec logo, A-dec Inspire, Cascade, Century Plus, Continental, Decade, ICX, ICV, Performer, Preference, Preference Collection, Preference ICC, Radius, and reliablecreativesolutions are trademarks of A-dec Inc. and are registered in the United States and other countries. A-dec 500, A-dec 400, A-dec 300, and A-dec 200 are also trademarks of A-dec Inc. None of the trademarks or trade names in this document may be reproduced, copied, or manipulated in any manner without the express, written approval of the trademark owner.

Certain touchpad symbols are proprietary to A-dec Inc. Any use of these symbols, in whole or in part, without the express written consent of A-dec Inc. is strictly prohibited.

The colors and products shown may vary from the actual materials and are subject to change without notice. Please see your authorized A-dec dealer for accurate color samples and the most current product information.



©2018 A-dec Inc. All rights reserved.
85.0093.00/PP/2.5M/10-18(10-18)/Rev E2